

## PROTOCOL

# Human IFN- $\gamma$ / Granzyme B Double-Color Enzymatic ELISPOT Assay

### CONTENTS

- Human IFN- $\gamma$  Capture Ab
- Anti-human IFN- $\gamma$  (FITC) Detection Ab
- FITC-HRP
- Human Granzyme B Capture Ab
- Anti-human Granzyme B (Biotin) Detection Ab
- Strep-AP
- CTL-Test™ Medium
- Diluent A
- Diluent B
- Diluent C
- Diluent Blue
- Diluent Red
- S1 (Blue substrate component 1)
- S2 (Blue substrate component 2)
- S3 (Blue substrate component 3)
- R1 (Red substrate component 1)
- R2 (Red substrate component 2)
- Plates: 96-well, high-protein-binding, PVDF filter plates
- Adhesive plate sealing sheet
- Protocol

### PROCEDURE

#### DAY 0 — STERILE CONDITIONS

- Prepare Human IFN- $\gamma$ /Granzyme B Capture Solution and prepare 70% ethanol (see Solutions).
- Remove plate underdrain, pipette 15 $\mu$ l of 70% ethanol into each well quickly. Add 150 $\mu$ l of PBS, decant, and wash with 150 $\mu$ l of PBS two more times. (If using strip plates, there is no underdrain to remove before prewetting.)  
**Note:** Activation of the membrane with ethanol is instantaneous and can be seen visually as a graying of the membrane. Ethanol should be washed off as quickly as possible following activation.
- Replace underdrain and immediately (before plate dries) pipette 80 $\mu$ l/well Human IFN- $\gamma$ /Granzyme B Capture Solution. Seal plate with parafilm and incubate at 4°C overnight.

#### DAY 1 — STERILE CONDITIONS

- Prepare CTL-Test™ Medium (see Solutions).
- Prepare antigen/mitogen solutions at two times final concentration in CTL-Test™ Medium.
- Decant plate containing Capture Solution from Day 0 and wash one time with 150 $\mu$ l PBS.
- Plate antigen/mitogen solutions, 100 $\mu$ l/well. Ensure the pH and temperature are ideal for cells by placing the plate containing antigens into a 37°C incubator, 5-9% CO<sub>2</sub> if for more than 10-20 minutes before plating cells.
- Adjust PBMC to desired concentration in CTL-Test™ Medium, e.g.: 3 million/ml corresponding to 300,000 cells/well (cell numbers can be adjusted according to expected spot counts since 100,000-800,000 cells/well will provide linear results). Keep cells at 37°C in humidified incubator, 5-9% CO<sub>2</sub> while processing PBMC and until plating.
- Plate PBMC, 100 $\mu$ l/well using large orifice tips. Once completed, gently tap the sides of the plate and immediately place into a 37°C humidified incubator, 5-9% CO<sub>2</sub>.
- Incubate for 24 hours. (Prestimulation of cells with an antigen is often needed to elicit a Granzyme B response. The duration of this prestimulation is antigen-dependent and can last up to one week. Alternatively, PHA is an excellent positive control and is only a 24-hour assay.) Do not stack plates. Avoid shaking plates by carefully opening and closing incubator door. Do not touch plates during incubation.

#### DAY 2

- Prepare Buffer Solutions: PBS, distilled water and Tween-PBS (see Wash Buffers).
- Prepare Anti-human IFN- $\gamma$ /Granzyme B Detection Solution (see Solutions).
- Wash plate two times with PBS and then two times with 0.05% Tween-PBS, 200 $\mu$ l/well each time.
- Add 80 $\mu$ l/well Anti-human IFN- $\gamma$ /Granzyme B Detection Solution. Incubate at room temperature, two hours.
- Prepare Tertiary Solution (see Solutions).
- Wash plate three times with 0.05% Tween-PBS, 200 $\mu$ l/well.
- Add 80 $\mu$ l/well of Tertiary Solution. Incubate at room temperature, one hour.
- During incubation, prepare Blue and Red Developer Solutions (see Solutions).
- Wash plate two times with 0.05% Tween-PBS, and then two times with distilled water, 200 $\mu$ l/well each time.
- Add Blue Developer Solution, 80 $\mu$ l/well. Incubate at room temperature, 15 minutes.
- Stop reaction by gently rinsing membrane with tap water, decant, and repeat three times.
- Decant plate with tap water and wash once more with distilled water, 200 $\mu$ l/well.
- Add Red Developer Solution, 80 $\mu$ l/well. Incubate at room temperature, 5-10 minutes.
- Stop reaction by gently rinsing membrane with tap water, decant, and repeat three times.
- Remove protective underdrain from the plate and rinse back of plate with tap water.
- Air-dry plate for two hours in running laminar flow hood or for 24 hours face down on paper towels on bench top.
- Scan and count plate. (CTL has scanning and analysis services available and offers a trial version of ImmunoSpot® Software with the purchase of any kit. Email kitscanningservices@immunospot.com for more info.)




Visit our YouTube channel for several helpful videos on working with ELISPOT assays and PBMC: [www.youtube.com/user/ImmunoSpot](http://www.youtube.com/user/ImmunoSpot).

The CTL Thawing Protocol for Cryopreserved Human PBMC is available at [www.immunospot.com](http://www.immunospot.com).

See other side for Solutions and Technical Tips.

For laboratory research use only. Not for use in diagnostic or therapeutic procedures.

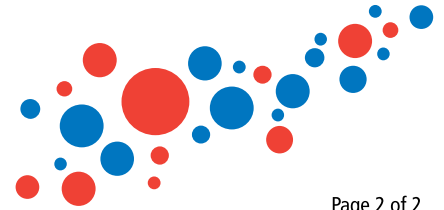

## SOLUTIONS

All solutions should be freshly-made prior to use. It is important to quick-spin the vials before use to ensure content volumes.

- **70% Ethanol (not included):** Dilute 190-200 proof ethanol. For 10ml, add 7ml of ethanol to 3ml of distilled water.
- **CTL-Test™ Medium:** Prepare medium by adding 1% fresh L-glutamine. The amount of medium needed will depend on variables such as cell yield and number of samples tested but will be no less than 20ml for one full plate; warm to 37°C before using.
- **Capture Solution:** Dilute *Human IFN- $\gamma$*  and *Human Granzyme B Capture Antibodies* in *Diluent A*. For one plate, add 40 $\mu$ l of *Human IFN- $\gamma$  Capture Antibody* and 80 $\mu$ l of *Human Granzyme B Capture Antibodies* to 10ml of *Diluent A*.
- **Detection Solution:** Dilute *Anti-human IFN- $\gamma$  (FITC)* and *Anti-human Granzyme B (Biotin) Detection Antibodies* in *Diluent B*. For one plate, add 40 $\mu$ l of *Anti-human IFN- $\gamma$  (FITC)* and 10 $\mu$ l of *Anti-human Granzyme B (Biotin) Detection Antibodies* to 10ml of *Diluent B* and filter through a 0.1 $\mu$ m low protein binding filter.
- **Tertiary Solution:** Dilute *FITC-HRP* 1:1000 and *Strep-AP* 1:1000 in *Diluent C*. For one plate, add 10 $\mu$ l of *FITC-HRP* and 10 $\mu$ l of *Strep-AP* to 10ml of *Diluent C*.
- **Blue Developer Solution:** To develop *Granzyme B* spots, add the *Substrate Solutions* in sequential steps to 10ml of *Diluent Blue*.  
For one plate:  
Step 1 – Add 160 $\mu$ l of *S1* to 10ml of *Diluent Blue*. Mix well!  
Step 2 – Add 160 $\mu$ l of *S2*. Mix well!  
Step 3 – Add 92 $\mu$ l of *S3*. Mix well!
- **Red Developer Solution:** To develop *IFN- $\gamma$*  spots, add the *Substrate Solutions* in sequential steps to 10ml of *Diluent Red*.  
For one plate:  
Step 1 – Add 180 $\mu$ l of *R1* to 10ml of *Diluent Red*. Mix well!  
Step 2 – Add 160 $\mu$ l of *R2*. Mix well!

It is recommended to make the Blue and Red Developer Solutions within ten minutes of use and keep them protected from direct light.

### Wash Buffers (not included)

#### For each plate prepare:

- 0.05% Tween-PBS: 100 $\mu$ l Tween-20 in 200ml PBS
- PBS, sterile, 100ml
- Distilled water, 100ml

## TECHNICAL TIPS

- Upon successful completion of the assay, IFN- $\gamma$  spots will be red, Granzyme B spots will be blue.
- CTL highly recommends doing single-color, positive control wells for color compensation during analysis.
- To maximize the use of each non-precoated plate, an adhesive plate-sealing sheet has been included that can be adhered to the top of the plate to cover unused wells for use in subsequent assays. Use your thumbs to firmly adhere the sheet to the plate and a razor blade to cut the sheet to expose only the necessary wells.
- We highly recommend the use of CTL Serum-free Media for freezing, washing, and testing PBMC. Even brief exposure to a mitogenic serum can cause high background while other sera can have suppressive effects. CTL also recommends using the CTL-LDC™ Kit for accurate live/dead cell counts.
- Deviations from specified temperatures, timing requirements, number of washing steps, and specified reagent preparation volumes may alter the performance of the assay.
- Low protein binding PVDF syringe filters (Millipore catalog #SLVV033RS) should be used for filtration of the Detection Solution to avoid loss of protein content.
- Plates may be washed manually or with a suitable automated plate washer with adjusted pin length and flow rate so membranes and spots are not damaged (CTL recommends the CTL 405LSR).
- To avoid damage to the PVDF membrane in the wells, do not touch the membrane with pipette tips or with the plate washer. The PVDF membrane is permeable and protected by an underdrain. Avoid direct contact between the well bottom and wet surfaces, including paper towels or any other materials that will absorb liquid.
- While processing plates, the PVDF membrane at the bottom of the wells must remain wet.
- After completion of the experiment, do not dry the ELISPOT assay plates at temperatures exceeding 37°C as this may cause the membrane to crack.
- Spots may not be readily visible while the membrane is still wet. Scan and count plates only after membranes have completely dried.
- Higher background appearing in the control wells can be potentially overcome using the following steps:
  - When working with precultured cells, wash the cells thoroughly in CTL-Wash™ prior to the experiment in order to avoid carryover of cytokines and other substances; use CTL-Test™ for testing PBMC.
  - The SmartCount™ module of the ImmunoSpot® counting software automatically recognizes spots over high background or uneven background, correcting background deviations. The Autogating™ module will help discern between T cell-derived and background spots. The CTL technical support team will gladly assist you with using the ImmunoSpot® Software for the analysis of complicated test results.
- Data analysis: The CTL ImmunoSpot® Analyzers along with the ImmunoSpot® Software have advanced features that permit automated, objective recognition of spots, gating and counting. An ELISPOT data management tool, SpotMap®, is also available to facilitate high-throughput ELISPOT work.

The CTL team will gladly assist you with data analysis and troubleshooting, as well as in customizing ELISPOT assays to suit your needs. Please contact us at [kits@immunospot.com](mailto:kits@immunospot.com).



Visit our YouTube channel for several helpful videos on working with ELISPOT assays and PBMC: [www.youtube.com/user/ImmunoSpot](http://www.youtube.com/user/ImmunoSpot).

The CTL Thawing Protocol for Cryopreserved Human PBMC is available at [www.immunospot.com](http://www.immunospot.com).

See other side for Contents and Procedure.

For laboratory research use only. Not for use in diagnostic or therapeutic procedures.



Cellular Technology Limited • 20521 Chagrin Blvd. • Shaker Hts., OH 44122-5350 USA  
+1 216-791-5084 Phone • +1 888-791-4005 Toll Free US • +1 216-751-1373 Fax  
[kits@immunospot.com](mailto:kits@immunospot.com) • [www.immunospot.com](http://www.immunospot.com)

CTL-Europe GmbH • Hans-Boeckler-Str. 19-29, 53225 • Bonn, Germany  
+49 (0) 228 8544 78-0 Phone • +49 (0) 228 8544 78-22 Fax  
[ctl-europe@immunospot.eu](mailto:ctl-europe@immunospot.eu) • [www.immunospot.eu](http://www.immunospot.eu)