




Product Data Sheet

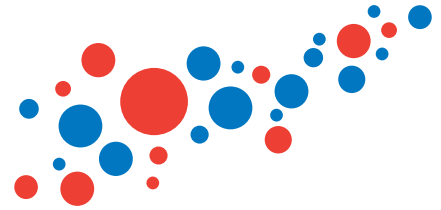

HCMV pp65 (495-503) (frozen, ready-to-use peptide solution)

Catalog #: CEF32-07-001

Product Name: HCMVpp65 (495-503)

Peptide Sequence: NLVPMVATV

Packaging Size: 1.5ml per vial, peptide dissolved in serum-free CTL-Test™ media

Concentration: 10nm

Description: HCMV pp65 (495-503) is an HLA-A2-restricted synthetic peptide derived from human cytomegalovirus.

Activity: HCMV pp65 (495-503) stimulates peptide-specific CD8+ memory T cells to produce IFN- γ .

Recommended Use: After thawing, the reagent is a 2x ready-to-use working solution (concentration 10nm). When working in 96 well plates, add 100 μ l of peptide solution to 100uL of cells per well. The cell concentration should be 100,000 – 800,000 cells/well.

Storage: Keep frozen at -20°C to -80°C.

References: Currier JR, Kuta EG, Turk E, Earhart LB, Loomis-Price L, Janetzki S, Ferrari G, Birx DL, Cox JH. A panel of MHC class I restricted viral peptides for use as a quality control for vaccine trial ELISPOT assays. *J Immunol Methods* 260:157-172, 2002.

Mwau M, Michael A and Hanke T. Design and validation of an enzyme-linked immunospot assay for use in clinical trials of candidate HIV vaccines. *AIDS Res Hum Retroviruses* 18: 611-8, 2002.

Warning: **For laboratory research use only, not for diagnostic use. Not for use in humans.**

Please note that the information contained in this document is privileged and confidential and protected from disclosure. You are hereby notified that any dissemination, distribution, or copying of this document is strictly prohibited. © 2012. All rights reserved.